

ASAHI Meister 16

ASAHI peripheral vascular guide wire

Torqueable + Shapeable

ASAHI Meister 16

ASAHI peripheral vascular guide wire

Torqueable

Precise torqueability
with ACT ONE
(One to one torque)

Shapeable

Stainless steel core allows
easy tip shaping and re-shaping
compared to nitinol core


Trackable

Polymer jacket and
well-designed shaft enhance
tracking ability


What is ACT ONE?



Employs the ASAHI brand rope coil,
which provides torqueability, flexibility and resiliency.

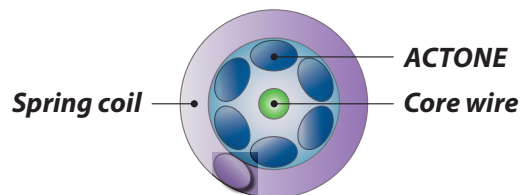


Enhanced torque response offers one-to-one torque control

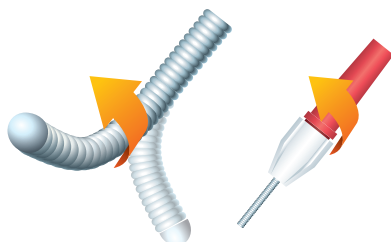

ASAHI Meister 16

Efficient transmission of proximal rotation to the distal tip.

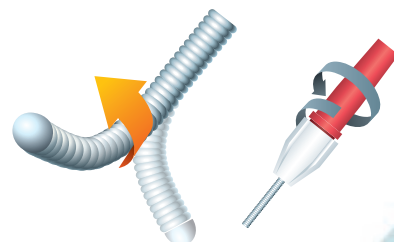


ASAHI conventional guide wire

Torque delay due to resistance causes whipping motion

peable + Trackable

0.41 mm (0.016 inch) guide wire
for peripheral vascular interventions

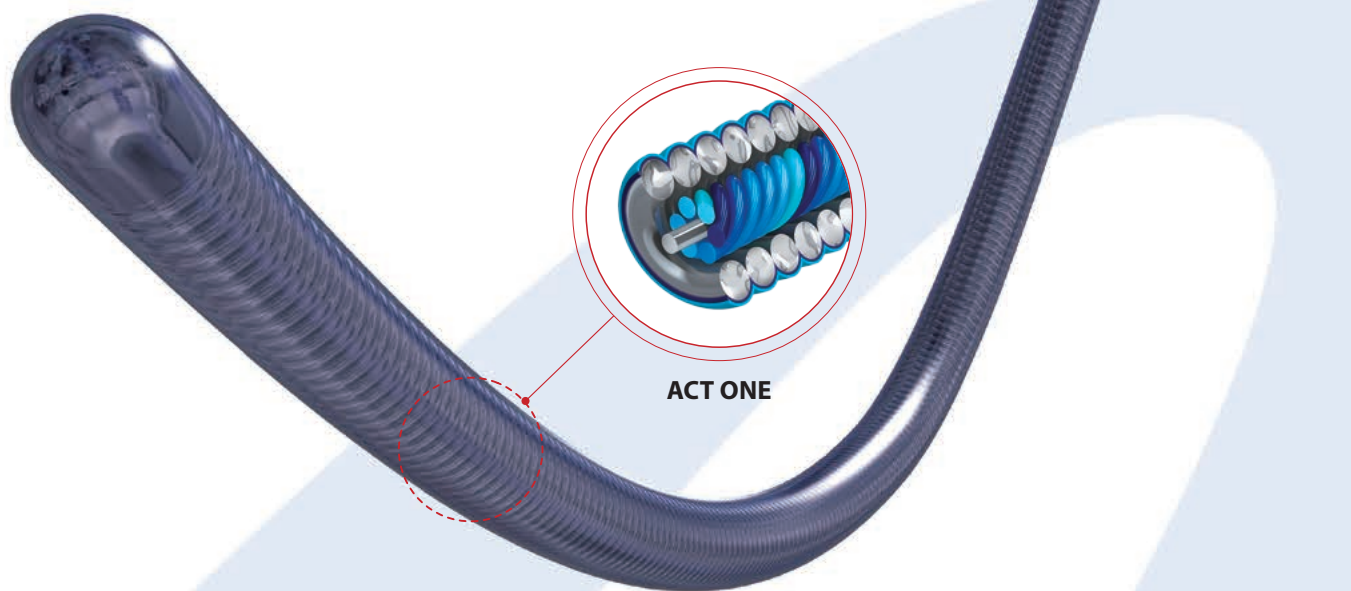


Pre-shaped tip

Thanks to pre-shaped ACT ONE, ASAHI Meister 16 has a good shape retention



Angled 45°

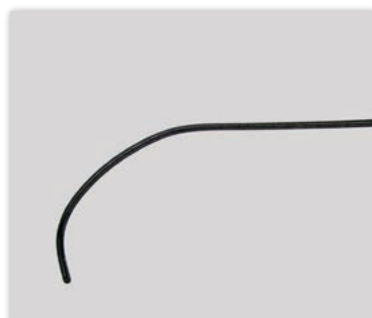


Angled double



Round curve

Structure & ordering information

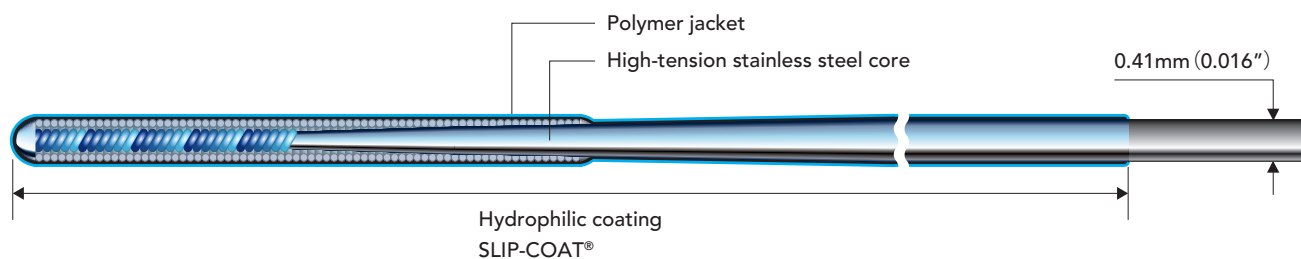

Product	Catalog No.	Usable length	Shaft O.D.	Radiopaque length	Polymer jacket length	Hydrophilic coating length	Tip shape
ASAHI Meister16	WAMS-165-1645	165 _{cm}	0.41 _{mm} (0.016 _{inch})	5 _{cm}	160 _{cm}	160 _{cm}	Angled 45°
	WAMS-180-1645	180 _{cm}	0.41 _{mm} (0.016 _{inch})	5 _{cm}	175 _{cm}	175 _{cm}	Angled 45°
	WAMS-165-16WA	165 _{cm}	0.41 _{mm} (0.016 _{inch})	5 _{cm}	160 _{cm}	160 _{cm}	Angled double
	WAMS-180-16WA	180 _{cm}	0.41 _{mm} (0.016 _{inch})	5 _{cm}	175 _{cm}	175 _{cm}	Angled double
	WAMS-165-16ST	165 _{cm}	0.41 _{mm} (0.016 _{inch})	5 _{cm}	160 _{cm}	160 _{cm}	Round curve
	WAMS-180-16ST	180 _{cm}	0.41 _{mm} (0.016 _{inch})	5 _{cm}	175 _{cm}	175 _{cm}	Round curve

Not all shapes and sizes are available in all markets. Please check availability with your sales representative.
Coated with SLIP-COAT® coating.

© 2018 ASAHI INTECC CO., LTD.
ASAHI, *Meister* and *ACT ONE* are trademarks or registered trademarks of ASAHI INTECC CO., LTD. in Japan and other countries.

CE 0344

Manufactured by:

Distributed by:

 **ASAHI INTECC**

ASAHI INTECC CO.,LTD.

3-100 Akatsuki-cho, Seto, Aichi 489-0071 Japan

Contact phone number: +81-561-86-9101

<http://www.asahi-intecc.com>