# Stiff CHIKAI black CHIKAI black 18 CHIKAI black 14 soft tip ASAHI CHIKAI black 18 soft tip ASAHI CHIKAI black 18 soft tip ASAHI CHIKAI black 18 soft tip Low Support High

# normal vessels / standard cases





tortuous vessels / complex cases

CHIK AI black

ASAHI
CHIK AI black 14 soft tip

# tiny vessels / precision cases





large vessels / high support cases


CHIK-AI black 18

CHIK-AI black 18

CHIK-AI black 18 soft tip

# Specification & ordering information

Product	Catalog No.	O.D.	Coating	Usable length	Coil length	Radiopaque <b>l</b> ength	Tip shape
ASAHI CHIKAI 008	WAIN-CKI-008-200	0.20/0.26 <sub>mm</sub> (0.008/0.010 <sub>inch</sub> )	Hydrophilic coating (SLIP-COAT®) 180 cm	200 cm	9 cm	9 cm	Straight
ASAHI CHIKAI 10	WAIN-CKI-10-200	0.26 <sub>mm</sub> (0.010 <sub>inch</sub> )	Hydrophilic coating (SLIP-COAT®) 170 cm	200 cm	9.5 cm	3 cm	Straight
	WAIN-CKI-10-300	0.26 mm (0.010 inch)	Hydrophilic coating (SLIP-COAT®) 170cm	300 cm	9.5 cm	3 cm	Straight
ASAHI CHIKAI	WAIN-CKI-200	0.36 <sub>mm</sub> (0.014 <sub>inch</sub> )	Hydrophilic coating (SLIP-COAT®) 170cm	200 cm	30 cm	5 cm	Straight
	WAIN-CKI-200-RC	0.36 mm (0.014 inch)	Hydrophilic coating (SLIP-COAT®) 170cm	200 cm	30 cm	5 cm	Round curve
	WAIN-CKI-300	0.36 <sub>mm</sub> (0.014 <sub>inch</sub> )	Hydrophilic coating (SLIP-COAT®) 170cm	300 cm	30 cm	5 cm	Straight
	WAIN-CKI-300-RC	0.36 mm (0.014 inch)	Hydrophilic coating (SLIP-COAT®) 170cm	300 cm	30 cm	5 cm	Round curve
ASAHI CHIKAI black	WAIN-CKI-200-BS	0.36 <sub>mm</sub> (0.014 <sub>inch</sub> )	Polymer jacket / Hydrophilic coating (SLIP-COAT®) 170 cm	200 cm	30 cm	<b>3</b> cm	Round curve
	WAIN-CKI-200-BA	0.36 <sub>mm</sub> (0.014 <sub>inch</sub> )	Polymer jacket / Hydrophilic coating (SLIP-COAT®) 150cm	200 cm	30 cm	3 cm	Angled 90°
ASAHI CHIKAI black 14 soft tip	AIN-CKI-200-B-SFT	0.36 <sub>mm</sub> (0.014 <sub>inch</sub> )	Polymer jacket / Hydrophilic coating (SLIP-COAT®) 170 cm	200 cm	30 cm	3 cm	Round curve
ASAHI CHIKAI black 18	WAIN-CKI-18-200-BS	0.45/0.36mm (0.018/0.014inch)	Polymer jacket / Hydrophilic coating (SLIP-COAT®) 170 cm	200 cm	34cm	<b>5</b> cm	Round curve
ASAHI CHIKAI black 18 soft tip	AIN-CKI-18-200-SFT	0.45/0.36mm (0.018/0.014inch)	Polymer jacket / Hydrophilic coating 150 cm / (SLIP-COAT®) 170 cm	200 cm	34cm	5 cm	Round curve

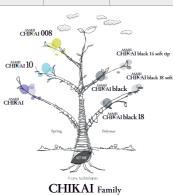

**( (** 0344

Not all shapes and sizes are available in all markets. Please check availability with your sales representative. Coated with  $SLIP-COAT^8$  coating.

CHIKAI" and "ACT ONE" are trademarks or registered trademarks of ASAHI INTECC CO., LTD. in Japa

Distributed by:

Manufactured by:



# ASAHI INTECC CO.,LTD.


3-100 Akatsuki-cho, Seto, Aichi 489-0071 Japan Contact phone number: +81-561-86-9101 http://www.asahi-intecc.com



# **ASAHI** HIKAI Family

# Providing good torque response, support and durability for neuro interventional procedures

ASAHI CHIKAI family of guide wires has been originally developed based on Asahi Intecc's unique set of technologies which accomplished great success in the field of cardiovascular interventions.

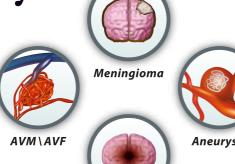





# **Technologies with Dual-coil design using ACT ONE**



Employs the ASAHI brand rope coil, which provides torqueability, flexibility and resiliency.

ASAHI INTECC CO., LTD.'s unique technologies in wire drawing, wire forming, torque transmission and coating are utilized to create highly responsive guide wires featuring ASAHI's ACT ONE, dual-coil design.



# **Key features**

- Precise torque response
- ACT ONE provides exceptional torque response in tortuous anatomy.
- **Durable tip-shape retention**
- Incorporating with ACT ONE, re-shapable pre-shaped round curve and 90 degree angle enhance tip shape memory and durability.








\*Available for CHIKAI CHIKAI black, CHIKAI black 18,

and CHIKAI black 14/18 soft tip

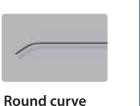








ASAHI CHIKAI, ASAHI CHIKAI 10 and ASAHI CHIKAI 008 are designed for an optimal balance between tip flexibility and microcatheter support

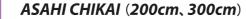

- 0.36mm (0.014") workhorse guide wire

- Tip radiopacity . SLIP-COAT® hydrophilic coating ..... 170 cm
  - ...5 cm Well-balanced workhorse wire



## ASAHI CHIKAI 008 (200cm)

Soft flexible tip & distal body

0.26mm (0.010")

....9 cm SLIP-COAT® hydrophilic coating ..... 180 cm

- Distal anatomy / small-sized vessels
- ▶ Compatible with microcatheters with small inner lumen





**ASAHI CHIKAI 10 (200cm, 300cm)** 

Balanced tip and higher body support

- ....3 cm Mid to distal anatomy / larger-sized vessels
- SLIP-COAT® hydrophilic coating ..... 170 cm ▶ Compatible with microcatheters with relatively larger inner lumen














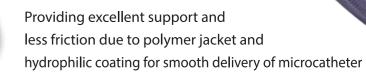

ASAHI CHIKAI black (200cm)

- 0.36mm (0.014") guide wire with enhanced lubricity

0.36mm (0.014") Tip radiopacity ... ...3 cm Complex cases and targeting distal anatomy

- Polymer jacket length ..... .150 cm ► Ease distal aspiration catheter placement
- SLIP-COAT® hydrophilic coating ..... 170 cm

ASAHI CHIKAI black 14 soft tip (200cm) Safe manipulation

- Tip radiopacity .... ..3 cm Complex cases and targeting distal anatomy Polymer jacket length ... .150 cm
- ▶ Ease distal aspiration catheter placement SLIP-COAT® hydrophilic coating ..... 170 cm ▶ Softer tip enhances safety in challenging anatomies
  - ▶ 43% softer tip than CHIKAI black
- **ASAHI CHIKAI black 18 (200cm)**

Speciality guide wire

0.41mm (0.016")

- ...5 cm Enough support to deliver large ID microcatheter Tip radiopacity ...
- Polymer jacket length ..... .. 150 cm Reverse tapered shaft design to avoid ledge effect ■ SLIP-COAT® hydrophilic coating ..... 170 cm ▶ Ease clot retriever delivery and positioning

ASAHI CHIKAI black 18 soft tip (200cm) Speciality guide wire

0.41mm (0.016"

...5 cm Enough support to deliver large ID microcatheter Polymer jacket length ... .. 150 cm

- ▶ Reverse tapered shaft design to avoid ledge effect ● SLIP-COAT® hydrophilic coating ..... 170 cm Ease clot retriever delivery and positioning
  - ▶ 30% softer tip than CHIKAI black 18