

ASAHI Gladius EX

Middleweight Non-Tapered Long Polymer Jacket Guide Wire

SPECIFICATIONS

ASAHI Gladius EX

TIP LOAD

3.0gf

COIL MATERIAL

Stainless Steel

CORE

Core + ACT ONE

WIRE O.D.

0.36mm (0.014 Inch)

COVER

Polymer Jacket

50cm Hydrophilic* COATING | 40cm Polymer Jacket COVER

3cm RADIOPACITY | 11cm SPRING COIL | Pre-Shaped TIP

ASAHI Gladius EX

ACT ONE

Catalogue No.	Diameter	Usable Length	Spring Coil Length	Hydrophilic* Coating Length	Radiopaque Length	Tip Shape
AP14R025P	0.36mm (0.014 Inch)	200cm	11cm	50cm	3cm	Pre-Shape
AP14R325P		300cm				

*Coated with SLIP-COAT® coating



© 2020 ASAHI INTECC CO., LTD.
"ASAHI", "ACT ONE" and "ASAHI Gladius" are trademarks or registered trademarks of ASAHI INTECC CO., LTD. in Japan and other countries.
Coated with SLIP-COAT® coating

Manufactured by:

ASAHI INTECC

ASAHI INTECC CO.,LTD.
3-100 Akatsuki-cho, Seto, Aichi 489-0071 Japan
Contact phone number: +81-561-86-9101
<http://www.asahi-intecc.com>


Distributed by:



ASAHI Gladius EX


ASAHI Gladius EX

Middleweight Non-Tapered Long Polymer Jacket Guide Wire

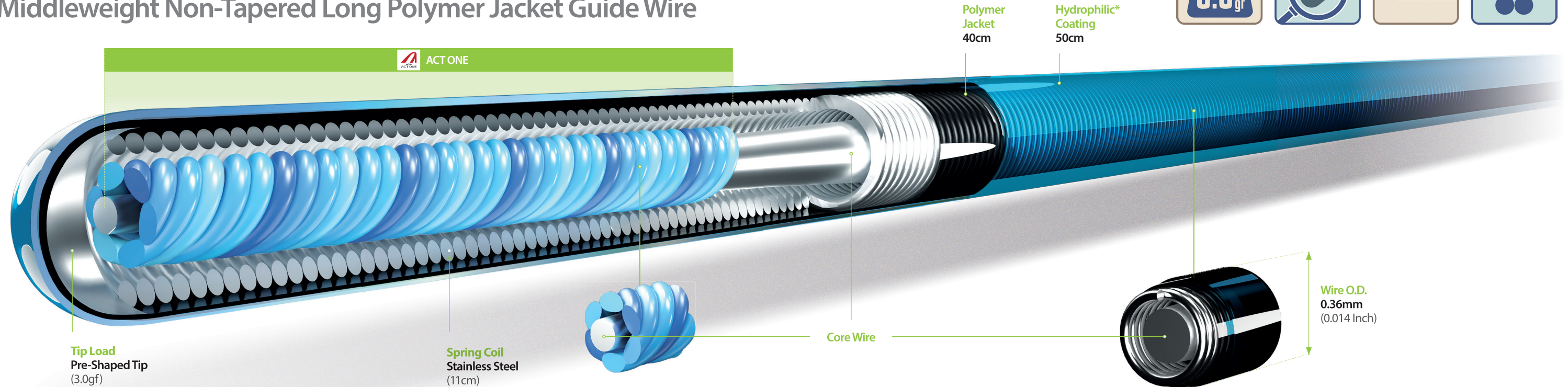

Mini Pre-Shape

1mm mini pre-shape manufactured to ensure tip shape within lesion which equals-easier tip shape control.

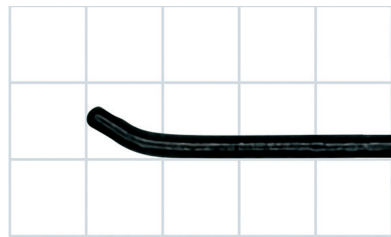

Long Coating Length

ASAHI Gladius EX has longer hydrophilic coating compared to competitor A to maintain its lubricity through long arteries.

ASAHI Gladius EX - 200cm




Competitor A



High Lubricity

ASAHI Gladius EX has excellent lubricity that allows the GW to pass smoothly through a long stenotic lesion

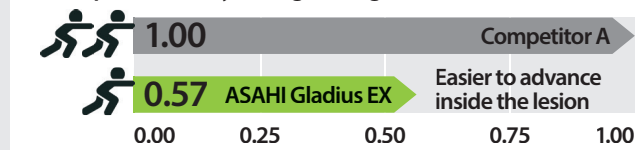


Lubricity Test Setup




1-To-1 Torque

ACT ONE

ASAHI Gladius EX has superior torque response without whip motion enabling excellent maneuverability

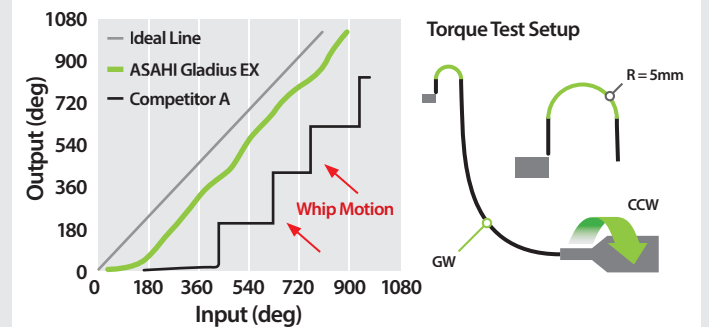

Tip Rigidity Comparison 1.5 - 30mm

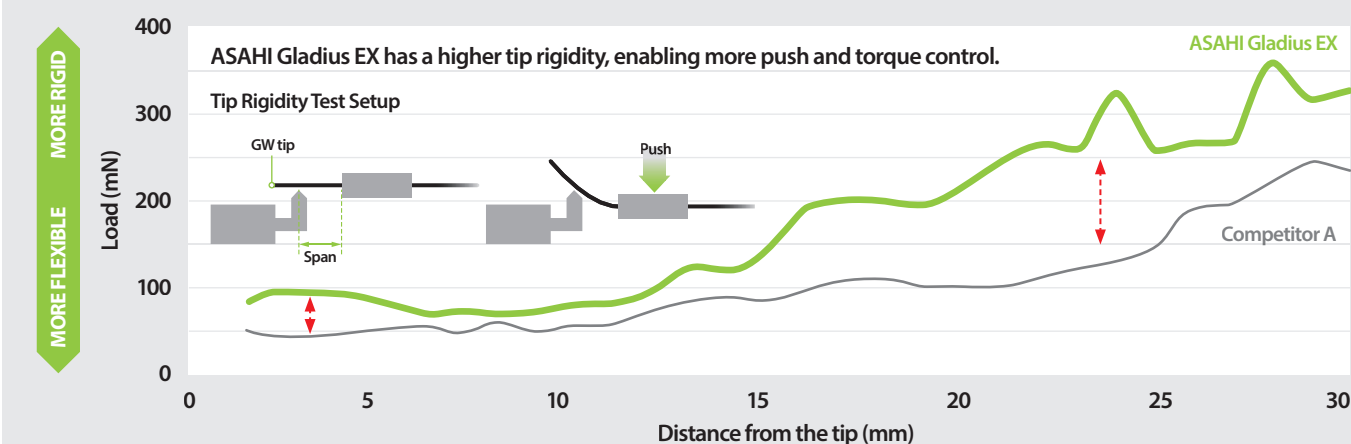

Shaft Rigidity Comparison 30 - 300mm

